




SAFETY STORAGE CABINETS

REMOVAL AND REPLACEMENT OF LOCKING MECHANISM (#25903)

- 1) REMOVE THE SCREW HOLDING THE ACCESS COVER PLATE. ROTATE THE PLATE TO GAIN ACCESS TO THE LOCKING MECHANISM (DO NOT REMOVE THE RIVET HOLDING THE COVER TO THE DOOR).
- 2) REMOVE THE RETAINER FROM THE SHAFT OF THE DOOR HANDLE CAREFULLY SO THAT LATCH RODS AND HARDWARE DO NOT FALL AWAY.
- 3) REMOVE CAM LEVER KEEPER, LATCH RODS, AND CAM LEVER/CAM LEVER LATCH FROM DOOR HANDLE SHAFT.
- 4) DOOR HANDLE CAN BE REMOVED BY UNSCREWING THE (2) SCREWS FROM THE INSIDE OF THE DOOR.
- 5) REPLACE WITH A NEW DOOR HANDLE AND (2) SCREWS.
- 6) SLIDE THE CAM LEVER/CAM LEVER LATCH COMPLETELY ONTO THE DOOR HANDLE SHAFT. CAREFULLY NOTE THE ORIENTATION OF THE PARTS USING THE DIAGRAM BELOW. THE CAM LEVER LATCH CAN HANG FREELY FOR THIS STEP.
- 7) TURN DOOR HANDLE 90 DEGREES TOWARDS THE HINGE.
- 8) GUIDE THE CAM LEVER LATCH INTO THE CAM GUIDE AND THROUGH THE SLOT IN THE DOOR.
- 9) REVERSE THE SEQUENCE OF DISASSEMBLY TO ASSEMBLE LATCH RODS, CAM LEVER KEEPER AND RETAINER.
- 10) ROTATE THE ACCESS COVER PLATE BACK INTO ITS ORIGINAL POSITION AND SECURE IT WITH THE SCREW.

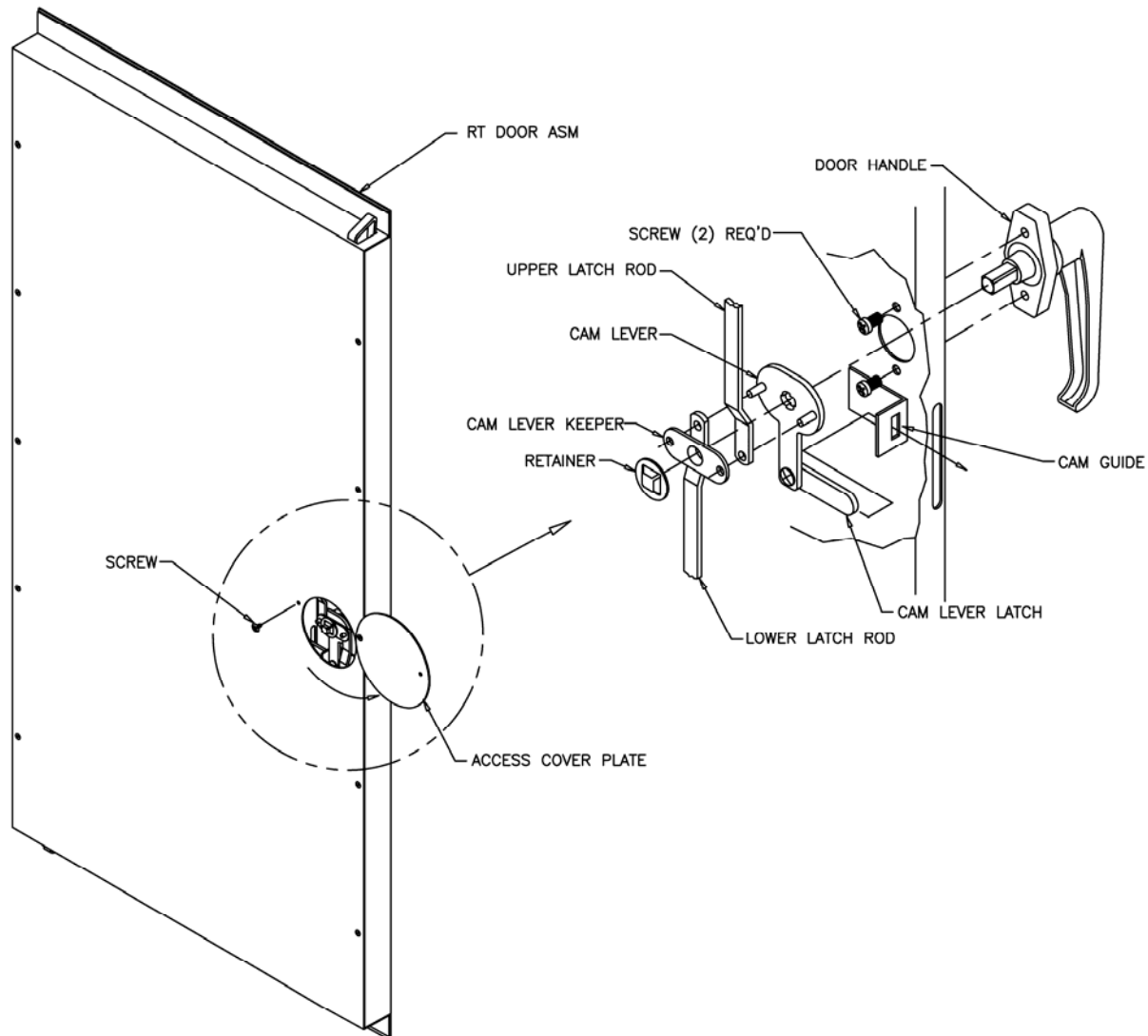

MANUFACTURER ASSUMES NO RESPONSIBILITY FOR INJURIES TO PERSONS OR DAMAGE TO CABINET, CONTENTS OR OTHER MATERIALS IF THE ABOVE INSTRUCTIONS ARE NOT COMPLIED WITH IN ALL PARTICULARS.